SET PÉČE O NOHY A PATY

10v1

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10in1


User manual

Please read the manual before first use. Keep the manual for future reference.

For hygienic reasons, grinding wheels are intended for one user only.

Description

The device is designed to remove hard skin from the feet. It leaves the feet silky soft and smooth. The foot treatment can be done both wet and dry. The device offers an ergonomic shape and design.

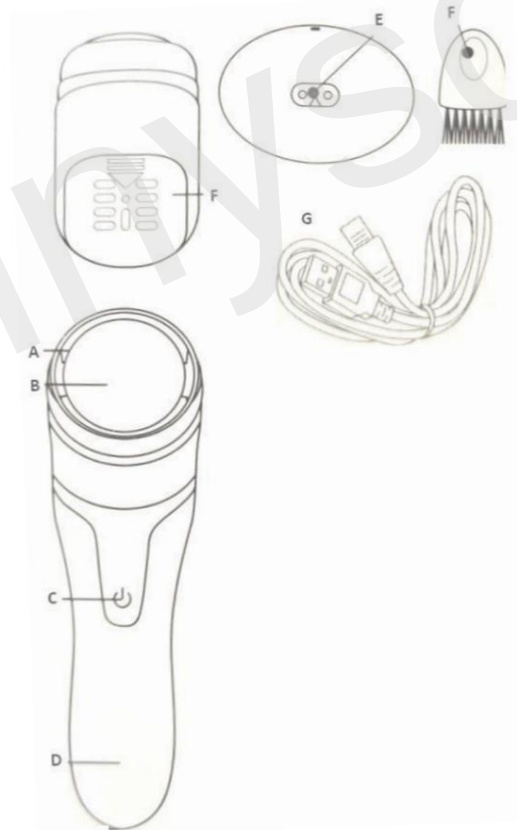- A – Protective cover
- B – Grinding head
- C – On/Off button + speed control
- D – LED indicator
- E – USB charging port
- F – Brush and dirt compartment
- G – USB charging cable

Warning/Important

1. Use the device only on the feet and heels.
2. The grinding wheel is intended for single use only (hygiene reasons). Do not share the wheels with other users.
3. Do not use the device if you suffer from diabetes or poor blood circulation.
4. If you suffer from a rash, eczema or other foot condition, consult your doctor before use.
5. Use the device on each spot for 2-3 seconds.
6. Do not use the device if the skin is red, pus-filled or blistered.
7. Do not use the device on areas where there are wounds, warts or ulcers.
8. If you experience pain during use, stop using the device immediately.
9. Excessive use may cause skin irritation.
10. If rash, irritation or bleeding persists, see a doctor.
11. The device is intended for indoor domestic use only. Do not use the device outdoors.
12. Keep the device out of reach of children. Keep the device out of reach of persons with mental, physical or psychological disabilities. Such persons may only use the device under the supervision of an instructed person.
13. Protect the device from direct sunlight.
14. Store the device in a dry and dust-free place at a temperature of 5-35 °C. Do not bend the charging cable.
15. Do not use the device without the grinding wheel attached.
16. Do not insert any objects into the device.
17. During use, protect your hair, nails, eyebrows, eyes, clothing, cables, shoelaces, etc. Otherwise, there is a risk of injury or damage to the device.
18. Do not press the area with too much force.
19. Do not disassemble or tamper with the device in any way.

Charging

Fully charge the device before using it for the first time. Charge the device for at least 2 hours. After a full charge, the device can be used for approximately 60 minutes. The time depends on the current conditions and skin condition. Charge, use and store the device at a temperature between 5-35 °C. Charge the device for approximately 2.5 hours.

Do not charge the device for more than 12 hours.

1. Turn off the device. Press the on/off button for 3 seconds.
2. Insert the charging cable into the device port. Make sure the cable is inserted in the correct orientation.

During charging, the battery status value will be displayed on the display. The device is fully charged when the value reaches 100%.

Use

1. First, cleanse and wash your skin.
2. Remove the protective cover of the grinding wheel. Make sure that the grinding head is correctly installed on the body of the device.
3. Press and hold the on/off button for 3 seconds to turn on the device. The device will turn on with gear 1 (ONE, default setting). Then press the on/off button briefly to switch to gear 2 (TWO).
4. Aim the head of the device at the desired location. Spend about 2-3 seconds on each location. Do not press the device too hard. Otherwise, the device will stop working. Stop the device after a while and check if the area is properly cleaned.
5. After finishing work, press and hold the on/off button for 3 seconds to turn off the device.
6. Then wipe the cleaned area with a damp towel.
7. Then clean the device head.
8. Finally, place the protective cover on the grinding head.

Note: Align the symbol on the cover  with the symbol  on the device.

Press and snap the cap into place.

Replacing the grinding head

1. Follow the picture below. Hold the holder and place the curved grooves under the grinding head. Then press the holder at both ends and pull. Remove the grinding head.
2. Then install the new head into the device. Insert it and push it in. You will hear a click. The head is then installed correctly. Make sure the head is installed correctly before use.

Maintenance and cleaning

Note: Clean the device after use. The device is waterproof. Remove the protective cover, sanding head and fold out the dirt compartment. Use a brush and then rinse it under running water. Do not turn on the device during cleaning.

Do not use salt water, water above 60°C, gasoline, alcohol and other corrosive liquids and detergents for cleaning. The device is waterproof with IPX6 protection.

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