BACtrack C8

Breathalyzer BT-C8 — User Guide

Select a language

EN	English.....	page 2
DE	Deutsch...	page 15
CS	Čeština ...	page 28
SK	Slovenčina	page 41
PL	Polski.....	page 54
HU	Magyar ...	page 68
RO	Română...	page 81
SL	Slovenščina	page 94
BG	Bulgarian	page 107

BACtrack C8

BT-C8

User Guide

Contents

1	Introduction	2
1.1	Alcohol and its effects on the human body	2
1.2	Alcohol is a drug	2
1.3	How does alcohol affect my body?	3
1.4	How does alcohol affect my driving?	3
2	Effects of alcohol*	4
3	Product overview	5
3.1	Important!	5
3.2	Components diagram	5
3.3	Operating modes	6
4	Standalone operation	6
5	BACtrack App operation	7
6	Troubleshooting	8
6.1	Alerts	8
6.2	FAQs	9
7	Precautions	9
7.1	Limits of breath alcohol testing	10
8	Specifications	10
9	Verification and calibration	10
9.1	Calibration instructions	11
10	Regulatory information	11
10.1	FCC regulatory compliance (USA)	11
10.2	RF exposure warning statement	11
10.3	Canada (ISED)	12
11	Service and support	12
11.1	Manufacturer	12

1. Introduction

This device is intended to measure alcohol in human breath. Measurements obtained by this device are used in the diagnosis of alcohol intoxication.

The BACtrack[®] is a professional alcohol screening device, used for the detection of alcohol in the breath. The BACtrack provides a digital result indicating the approximate BAC (Blood Alcohol Content) of the test subject. The BACtrack is easy to use: simply turn on the unit and begin testing in a matter of seconds.

Before you begin testing, please read this manual in its entirety and familiarize yourself with the BACtrack.

1.1. Alcohol and its effects on the human body

Alcohol is absorbed from the mouth, throat, stomach and intestines into the bloodstream.

Alcohol ingested by the human body can be detected in the breath because of its presence in the bloodstream.

Alcohol cannot be digested and it cannot be chemically changed in the bloodstream. As the blood flows through the lungs, alcohol in the bloodstream moves across the membranes of the lungs' air sacs (alveoli) into the air. The concentration of the alcohol in the alveolar air is directly related to the concentration of the alcohol in the blood. As the alveolar air is exhaled, the alcohol in it can be detected by the breath alcohol testing device.

The alcohol concentration in the breath is related to that in the blood, and because of this, an individual's BAC can be determined by measuring alcohol in the breath. The ratio of breath alcohol to blood alcohol is generally estimated to be 2,100:1. Therefore, 2,100 milliliters (ml) of alveolar air will contain approximately the same amount of alcohol as 1 ml of blood.

1.2. Alcohol is a drug

Alcohol is the chemical (ethanol or ethyl alcohol) resulting from the fermentation of grapes or grain. Alcohol is absorbed directly into your bloodstream. From your bloodstream, the fluids in your body tissues absorb the alcohol. Your brain is made up of a large concentration of fluids and will absorb a large amount of any alcohol you drink. Your liver eliminates the absorbed alcohol at its own rate and nothing you eat or drink can speed up the process. Your liver does this by oxidizing the alcohol (converting it into water and carbon dioxide). Coffee, food or any other "cure" will not sober you up; only time will do that.

1.3. How does alcohol affect my body?

Alcohol is a depressant. It has a relaxing effect on the muscles in your body. The muscles of your eyes relax and lose focus. Your eyesight will become fuzzy and you may experience double vision. Brain activity is slowed. Your judgment, reflexes and coordination are all negatively affected.

Some vision impairments that occur when you have been drinking include:

- narrowing of your field of vision;
- reduction in your depth perception;
- decreased ability to see in darkness;
- increased sensitivity to glare and a longer time for your eyes to readjust from the glare.

Some mental impairment can occur when you have been drinking, including:

- reduced awareness of danger;
- becoming overly confident and reckless;
- difficulty in making decisions;
- reduction in balance;
- slowed reflexes;
- impaired judgment.

1.4. How does alcohol affect my driving?

We can't emphasize enough how dangerous it is for you to drink and drive. Your vision and your brain are the most important factors in driving safely. If you jeopardize either by drinking, and then you drive, you are likely to be involved in a serious or fatal collision.

Do not use your BACtrack as a tool to determine whether you should operate a motor vehicle or equipment, or perform any other dangerous act. Do not drink and drive. Always have a designated driver when alcohol is being consumed.

Source: California Department of Motor Vehicles (<http://www.dmv.ca.gov>).

2. Effects of alcohol*

BAC	Dose-specific effects
0.02–0.03 %BAC	Slight euphoria. Loss of shyness. Depressant effects are not apparent. Impairment possible in some individuals.
0.04–0.06 %BAC	Feeling of well-being, relaxation, lower inhibitions and sensation of warmth. Euphoria. Some minor impairment of reasoning and memory. Lowering of caution. Driving skills may be impaired at this level of intoxication.
0.07–0.09 %BAC	Slight impairment of balance, speech, vision, reaction time and hearing. Euphoria. Judgment and self-control are reduced. Caution, reason and memory are impaired. Driving skills are always impaired at this level of intoxication and higher.
0.10–0.12 %BAC	Significant impairment of motor coordination and loss of good judgment. Speech may be slurred. Balance, vision, reaction time and hearing will be impaired. Euphoria.
0.13–0.15 %BAC	Gross motor impairment and lack of physical control. Blurred vision and major loss of balance. Euphoria is reduced and dysphoria (anxiety, restlessness) begins to appear.
0.16–0.20 %BAC	Dysphoria predominates. Nausea may appear.
0.25 %BAC	Need for assistance in walking. Total mental confusion. Dysphoria with nausea and some vomiting.
0.30 %BAC	Loss of consciousness. Onset of coma. Possible death due to respiratory arrest**.

The generally accepted legal standard for alcohol intoxication in the United States is 0.08 %. However, your driving skills can be impaired at any level above 0.00 %BAC. It is never safe to drink any amount of alcohol and drive.

* The effects of alcohol intoxication are greatly influenced by individual variations among users. Some users will be intoxicated at a much lower BAC than shown in the table above.

** Death can occur at a lower BAC in some individuals.


3. Product overview

3.1. Important!

For most accurate results, it's important to **wait at least 20 minutes after smoking, eating or drinking before testing**. This is a universal principle for all breathalyzers, not just yours. The reason is that if you have recently consumed alcohol, the residual alcohol in your mouth is very high and can cause a false reading.

3.2. Components diagram

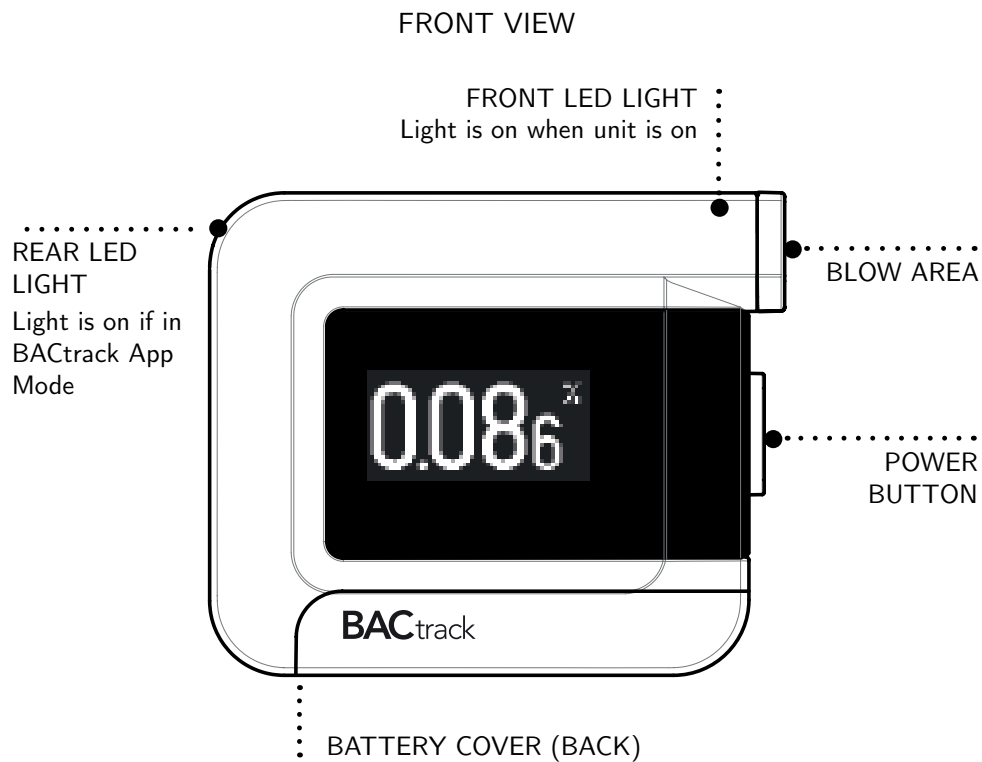3.3. Operating modes

You control how your BAC results are delivered, either:

1. to the screen directly on your BACtrack – this is **Standalone mode**. No results are saved or transmitted in this mode; or
2. to the BACtrack App on your mobile device – this is **BACtrack App mode**. You have the option to save results in this mode.

Note: The blue LED light on the rear of your BACtrack indicates which mode you are in. If the light is on, you are in BACtrack App mode. If the light is off, you are in Standalone mode. To switch from one mode to the other, hold the power button down for seven seconds.

4. Standalone operation

- STEP 1:** Install one AAA alkaline battery in the battery compartment.
- STEP 2:** If you would like to use a mouthpiece, attach it to the breathalyzer. Note: using a mouthpiece is optional.
- STEP 3:** Turn on your BACtrack by pressing the power button for one second. The front blue LED light indicates when it is on.
- STEP 4:** The BACtrack will display **Warming Up** and then count down. When the countdown approaches 2, begin to inhale a deep breath.
- STEP 5:** When the countdown reaches zero, **BLOW** will flash on the screen. Blow steadily and consistently for six seconds through the mouthpiece – do not blow as hard as possible. The BACtrack emits a beeping sound during the breath sample. Stop blowing when the beeping sound stops. There is a click at the end – that’s the air pump pulling air across the sensor.
- STEP 6:** Once the sensor has analyzed the breath sample, the estimated BAC value is displayed. Next, the ZeroLine feature will appear, showing approximately how long it will take for your estimated BAC level to return to 0.00 % if no additional drinks are consumed. The estimated time is displayed in hours and minutes. The screen will then display a brief summary of how your alcohol level may be affecting your body. The BACtrack will turn off automatically in 10 seconds. You can also turn off the breathalyzer during the countdown by pressing and holding down the power button for 2 seconds. The front blue LED light will turn off.

5. BACtrack App operation

- STEP 1:** Install one AAA alkaline battery in the battery compartment.
- STEP 2:** If you would like to use a mouthpiece, attach it to the breathalyzer.
Note: using a mouthpiece is optional.
- STEP 3:** Download the BACtrack App from the App Store or Google Play Store.
- STEP 4:** Go to Settings on your mobile device to confirm Bluetooth® is on; this will ensure your BACtrack will pair with the app.
- STEP 5:** Open the BACtrack App on your mobile device.
- STEP 6:** Press and hold the power button on your BACtrack for one second to turn on the device. The front blue LED light will turn on. If the rear blue LED light is off, hold the power button for seven seconds to switch to BACtrack App mode. Once in BACtrack App mode, the app and device will pair automatically. When the app pairs with the BACtrack, the app screen will momentarily show an image of your BACtrack.
- STEP 7:** Tap **BAC Test** to initiate the testing process. While the BACtrack sensor warms up, you have the opportunity to **Guess Your BAC**. This is optional and will not affect your final reading. Your BAC guess and actual result will be displayed so you can compare your perception to an actual number. Using this feature over time, you'll notice you are much better at understanding how alcohol affects your body.
- STEP 8:** When the screen on your mobile device reads **Blow Now**, take a deep breath and blow into your BACtrack continuously until all the blue segments are orange.
- STEP 9:** Once your BACtrack has analyzed the breath sample, it will transmit the reading to your mobile device, where the BAC value will be displayed.

STEP 10: Press **OK** to view the ZeroLine® feature and other information about your BAC result. ZeroLine will tell you the approximate time when your estimated BAC level will return to 0.00 % if no additional drinks are consumed.

Your BACtrack will power off automatically in 10 seconds. You can also turn it off at any time by pressing the power button for two seconds.

Please keep your BACtrack App updated as we continuously add new features and options.

6. Troubleshooting

6.1. Alerts

LOW BATTERY

Install a new alkaline AAA battery.

ERROR: TEMPERATURE

If the unit is powered on outside of the acceptable temperature range (0–40 °C or 32–104 °F), testing cannot be performed. Power off the unit. Retest once the BACtrack has been stored in an acceptable environment.

ERROR: BLOW STEADILY FOR 6 SECONDS

Please press the power button to restart the countdown cycle and retest.

ERROR: NO BLOW DETECTED

Please press the power button to restart the countdown cycle and retest.

ERROR: MAX BAC EXCEEDED

The estimated BAC exceeds the maximum measurable level of 0.30 %. If you did not wait 20 minutes after your last drink, any residual alcohol in your mouth could create a false high reading.

SEARCHING FOR DEVICE

The BACtrack is in BACtrack App mode and is searching for the BACtrack App on a mobile device. You can toggle between Standalone and BACtrack App mode by pressing the power button for 7 seconds.

NO DEVICE FOUND

The BACtrack is in BACtrack App mode and cannot locate the BACtrack App on a mobile device. Please make sure Bluetooth is enabled on the mobile device.

6.2. FAQs

Why is the rear blue LED light on?

You are in BACtrack App mode. Your results will be delivered to your BACtrack App.

Why are my results not showing up on the LED screen?

You are not in Standalone mode. Hold the power button for seven seconds to switch modes.

Why can't the app find my device?

Your BACtrack is not in BACtrack App mode. Hold the power button for seven seconds to switch modes.

7. Precautions

1. Wait at least 20 minutes after drinking, eating or smoking before testing. Not observing this waiting period can cause inaccurate readings and damage the sensor.
2. Do not blow smoke, food or liquids into the BACtrack, because this will damage the sensor.
3. Do not test in areas with strong winds, smoke, or in an area where a large amount of alcohol is being consumed.
4. Avoid testing in the presence of any substances that contain methyl alcohol, isopropyl alcohol or acetone. These substances may interfere with the results of the test.
5. Replace the AAA alkaline battery when the battery indicator icon reaches one bar. **Caution:** risk of explosion if the battery is replaced by an incorrect type. Dispose of the used battery according to the instructions.
6. Send your BACtrack in for periodic calibration service as required (see section 9).
7. Your BACtrack is designed to be used in a temperature range of 32–104 °F (0–40 °C).
8. The BAC results from this product are not admissible in a court of law.
9. Do not use your BACtrack as a tool to determine whether you should operate a motor vehicle or equipment, or perform any other dangerous act.
10. Do not drink and drive. Always have a designated driver when alcohol is being consumed.

7.1. Limits of breath alcohol testing

- This product cannot determine whether, or to what extent, a test subject's mental or physical functioning or judgment is impaired by alcohol, including whether the test subject is intoxicated under any definition of that word.
- Any act or failure to act based on a reading from this product is at the test subject's or user's own risk.
- The correlation between breath alcohol content and blood alcohol concentration depends on many variables, including environmental factors (such as air quality, wind, humidity and temperature) and the health condition of the test subject.
- A low BAC reading does not mean that the test subject's physical or mental performance or judgment (including reaction times) is sufficient to respond to an emergency.
- The concentration of alcohol in the blood of a test subject cannot be exactly determined by using a breath alcohol screening product.

8. Specifications

Dimensions	2.21 × 1.875 × 0.66 inches (5.62 × 3.45 × 1.68 cm)
Weight	2 oz (56 g), with mouthpiece and battery
Detection range	0.000–0.300 %BAC
Power supply	One AAA alkaline battery
Battery life	Approximately 300 tests
Warm-up time	10–20 seconds
Response time	10–20 seconds
Operating temperature	32–104 °F (0–40 °C)

9. Verification and calibration

The BACtrack is calibrated during manufacture using advanced alcohol simulation equipment. Known alcohol concentrations are passed through the sensor to set baseline values for testing. During testing, the unit compares breath samples to these baseline values.

The accuracy of breath alcohol testers can fluctuate after six to twelve months of normal use, depending on the number of tests performed and operating conditions.

If this product is providing inconsistent test results, is not providing any test results, or provides unusually high or low test results, the product should be recalibrated immediately by a trained professional.

The product should be recalibrated at least every six to twelve months. If the product is used every day, it may need to be recalibrated as frequently as every month.

9.1. Calibration instructions

Recalibration is arranged through the store where you bought your BACtrack or through the distributor (see section 11). Information on how to get your BACtrack recalibrated is also available at www.bactrack.com/calibration or, in the USA, by phone at 877.334.6876.

During recalibration, a professional technician will reset the sensor baseline values and also perform a full diagnostic check to ensure the accuracy and reliability of the product.

10. Regulatory information

10.1. FCC regulatory compliance (USA)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

The FCC requires the user to be notified that any changes or modifications made to this device that are not expressly approved by KHN Solutions Inc. may void the user's authority to operate the equipment.

10.2. RF exposure warning statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. End users must follow the specific operating instructions for satisfying RF exposure compliance.

10.3. Canada (ISED)

RSS-Gen, Issue 4, November 2014: This device complies with ISED licence-exempt RSS standard(s). Operation is subject to the following two conditions:

1. this device may not cause harmful interference;
2. this device must accept any interference received, including interference that may cause undesired operation of the device.

Radiation exposure statement: This equipment complies with the IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

French version of the ISED statements, as printed in the original:

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The information in this guide may change without notice. The manufacturer assumes no responsibility for any errors that may appear in this guide.

11. Service and support

If you believe your BACtrack is not operating properly, please first read this manual carefully. For further help, contact the store where you bought the product or the distributor (see the last page), or visit www.bactrack.com. Repairs and calibration are arranged through the distributor.

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