BACtrack S80

Breathalyzer — User Guide

Select a language

EN	English.....	page 2
DE	Deutsch...	page 13
CS	Čeština ...	page 24
SK	Slovenčina	page 34
PL	Polski.....	page 45
HU	Magyar ...	page 56
RO	Română...	page 67
SL	Slovenščina	page 78
BG	Bulgarian .	page 89

BACtrack® S80

User Guide

Contents

1	Introduction	1
1.1	Alcohol and Its Effects on the Human Body	1
1.2	Alcohol Is a Drug	1
1.3	How Does Alcohol Affect My Body?	2
1.4	How Does Alcohol Affect My Driving?	2
2	Effects of Alcohol*	3
3	Preparation	3
3.1	Installing the Batteries	3
3.2	Components Diagram	4
4	Operation	4
4.1	Error Notifications	5
4.2	Advanced Menu Options	6
5	Specifications	7
6	Precautions	7
7	Calibration	8
7.1	Calibration Instructions	8
8	Service and Support	8
8.1	Manufacturer	9

1. Introduction

This device is intended to measure alcohol in human breath. Measurements obtained by this device are used in the diagnosis of alcohol intoxication.

The BACtrack[®] is an alcohol screening device used for the detection of alcohol in the breath. The BACtrack provides a digital result indicating the approximate BAC (Blood Alcohol Content) of the test subject. The BACtrack is easy to use: simply turn on the unit and begin testing in a matter of seconds.

Before you begin testing, please read this manual in its entirety and familiarize yourself with the BACtrack.

1.1. Alcohol and Its Effects on the Human Body

Alcohol is absorbed from the mouth, throat, stomach and intestines into the bloodstream.

Alcohol ingested by the human body can be detected in the breath because of its presence in the bloodstream. Alcohol cannot be digested and it cannot be chemically changed in the bloodstream. As the blood flows through the lungs, alcohol in the bloodstream moves across the membranes of the lung's air sacs (alveoli) into the air. The concentration of the alcohol in the alveolar air is directly related to the concentration of the alcohol in the blood. As the alveolar air is exhaled, the alcohol in it can be detected by the breath alcohol testing device.

The alcohol concentration in the breath is related to that in the blood, and because of this, an individual's BAC can be determined by measuring alcohol in the breath. The ratio of breath alcohol to blood alcohol is generally estimated to be 2,100:1. Therefore, 2,100 milliliters (ml) of alveolar air will contain approximately the same amount of alcohol as 1 ml of blood.

1.2. Alcohol Is a Drug

Alcohol is the chemical (ethanol or ethyl alcohol) resulting from the fermentation of grapes or grain. Alcohol is absorbed directly into your bloodstream. From your bloodstream, the fluids in your body tissues absorb the alcohol. Your brain is made up of a large concentration of fluids and will absorb a large amount of any alcohol you drink. Your liver eliminates the absorbed alcohol at its own rate and nothing you eat or drink can speed up the process. Your liver does this by oxidizing the alcohol (converting it into water and carbon dioxide). Coffee, food or any other "cure" will not sober you up; only time will do that.

1.3. How Does Alcohol Affect My Body?

Alcohol is a depressant. It has a relaxing effect on the muscles in your body. The muscles of your eyes relax and lose focus. Your eyesight will become fuzzy and you may experience double vision. Brain activity is slowed. Your judgment, reflexes and coordination are all negatively affected.

Some vision impairments that occur when you have been drinking include:

- narrowing of your field of vision;
- reduction in your depth perception;
- decreased ability to see in darkness;

- increased sensitivity to glare and a longer time for your eyes to readjust from the glare.

Some mental impairment can occur when you have been drinking, including:

- reduced awareness of danger;
- becoming overly confident and reckless;
- difficulty in making decisions;
- reduction in balance;
- slowed reflexes;
- impaired judgment.

1.4. How Does Alcohol Affect My Driving?

We can't emphasize enough how dangerous it is for you to drink and drive. Your vision and your brain are the most important factors in driving safely. If you jeopardize either by drinking, and then you drive, you are likely to be involved in a serious or fatal collision.

Do not use your BACtrack as a tool to determine whether you should operate a motor vehicle or equipment, or perform any other dangerous act. Do not drink and drive. Always have a designated driver when alcohol is being consumed.

Source: California Department of Motor Vehicles, USA (<http://www.dmv.ca.gov>)

2. Effects of Alcohol*

%BAC	Effects
0.02–0.03%	Slight euphoria. Loss of shyness. Depressant effects are not apparent. Impairment possible in some individuals.
0.04–0.06%	Feeling of well-being, relaxation, lower inhibitions and sensation of warmth. Euphoria. Some minor impairment of reasoning and memory. Lowering of caution. Driving skills may be impaired at this level of intoxication.
0.07–0.09%	Slight impairment of balance, speech, vision, reaction time and hearing. Euphoria. Judgment and self-control are reduced. Caution, reason and memory are impaired. Driving skills are always impaired at this level of intoxication and higher.
0.10–0.12%	Significant impairment of motor coordination and loss of good judgment. Speech may be slurred. Balance, vision, reaction time and hearing will be impaired. Euphoria.
0.13–0.15%	Gross motor impairment and lack of physical control. Blurred vision and major loss of balance. Euphoria is reduced and dysphoria (anxiety, restlessness) begins to appear.
0.16–0.20%	Dysphoria predominates. Nausea may appear.
0.25%	Need for assistance in walking. Total mental confusion. Dysphoria with nausea and some vomiting.
0.30%	Loss of consciousness. Onset of coma. Possible death due to respiratory arrest**.

* The effects of alcohol intoxication are greatly influenced by individual variations among users. Some users will be intoxicated at a much lower BAC than shown above.

** Death can occur at a lower BAC in some individuals.


In Canada and the United States, the generally accepted legal standard for alcohol intoxication is 0.08% BAC; other countries set their own limits. However, your driving skills can be impaired at any level above 0.00% BAC. **It is never safe to drink any amount of alcohol and drive.**

3. Preparation

3.1. Installing the Batteries

Install two AA alkaline batteries in the battery compartment.

3.2. Components Diagram

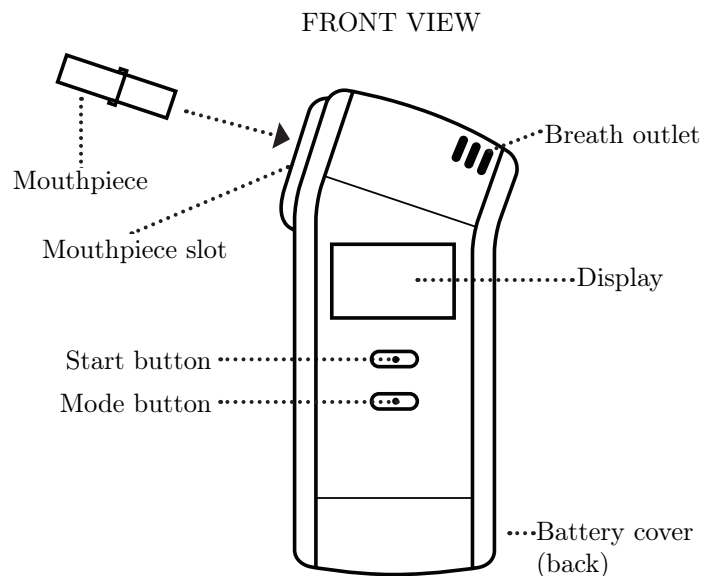
4. Operation


Test subjects should wait 20 minutes after eating, drinking or smoking before blowing into the BACtrack breathalyzer. Failure to wait 20 minutes can provide inaccurate test results and damage the sensor.

STEP 1: Insert a mouthpiece into the mouthpiece slot. Do not force it: the mouthpiece will only insert a few millimeters.

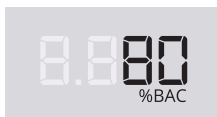STEP 2: Press **START**. The BACtrack will quickly display the total number of tests performed and then begin a countdown.

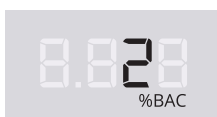STEP 3: When the countdown approaches 2, begin to take a deep breath.


STEP 4: When the countdown reaches zero, the **START** indicator flashes and the screen shows three lines. Blow through the mouthpiece for five seconds, until you hear a double beep.


STEP 5: Once the sensor has analyzed the breath sample, the estimated BAC value flashes for about 10 seconds.

The BACtrack breathalyzer turns off automatically after 3 minutes. You can also turn it off during the countdown by holding down **START** for 3 seconds.

4.1. Error Notifications

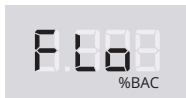
If the user does not blow a sufficient breath sample, the display shows FLo. Press **START** to restart the countdown cycle and test again.

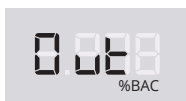
If no breath sample is blown within 3 minutes, the display shows Out. Press **START** to restart the countdown cycle and test again.


If the unit is switched on outside the acceptable temperature range (0–40 °C or 32–104 °F), it displays OutF and testing cannot be performed.

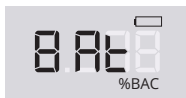
If the battery indicator is empty, install two new AA alkaline batteries.

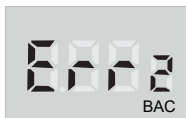If the unit displays Err2 after step 4, excessive saliva went into the sensor during the test. Wait 20 minutes, then try again.

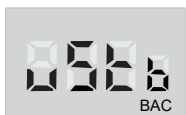This error indicates there may be environmental alcohol present which could interfere with the accuracy of the test results. Make sure there are no other sources of ethanol present (such as lotions or perfumes), wait at least one hour, then try again.

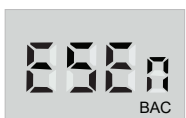This error indicates the sensor is not detecting normal at-rest values. Make sure your BACtrack has been stored in appropriate conditions, wait 20 minutes, then try again. If this error occurs consistently, contact support@bactrack.com.


This error indicates there may be physical sensor damage. Note that smoke or saliva can damage the sensor. Wait 10 minutes, then try again. If this error occurs consistently, contact support@bactrack.com.

4.2. Advanced Menu Options


The BACtrack offers several user-adjustable features through Menu Mode. To enter Menu Mode, turn on the BACtrack and, during the countdown, hold down **MODE** for 5 seconds. To move between the Menu Mode options, press **MODE**.


Units of measure: By default, the BACtrack displays test results in %BAC, which is the standard in North America and in many other regions. To change to a different unit of measure, press **START**. You can select per mille (‰) or mg/L as alternative units of measure.

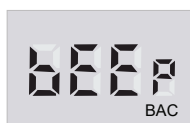Audible warnings: By default, the BACtrack does not sound an audible warning after a specific %BAC result. To set an audible warning after a specific %BAC test result, press **START** to cycle through the %BAC values.

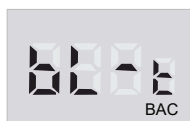
Blow time: The default blow time is 5 seconds. This longer blowing time generally provides a deep-lung air sample and a more accurate test result. To adjust the blow time, press **START**.

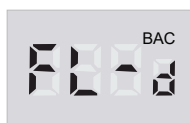Flow detection settings: Use this option to adjust how hard the user is required to blow. Press **START** to change the setting:


Standard flow detection level: This is the default setting of the unit. This level requires a less forceful breath from the user.

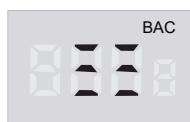Medium flow detection level: Press **START** to change to the medium flow level.

High flow detection level: Press **START** to change to the high flow level. This level requires a more forceful breath from the user.

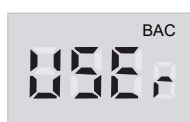User-activated test mode: Press **START** again to change to the user-activated test mode (USER). When testing in this mode, the operator must press **START** during a breath sample to activate the pump and perform a test.

5. Specifications

Dimensions	2.3 × 4.8 × 0.8 inches (5.8 × 12.2 × 2.2 cm)
Weight	4.8 oz (136 g) with mouthpiece and batteries
Sensor technology	Xtend™ fuel cell sensor technology
Detection range	0.000–0.500 %BAC
Test count	Displays the total number of tests performed
Power supply	Two AA alkaline batteries, included
Battery life	Approximately 1,000 tests
Warm-up time	10–30 seconds, depending on last use
Response time	3 seconds
Sensor accuracy	±0.005 %BAC at 0.050 %BAC
Operating temperature	32–104 °F (0–40 °C)

This device has been tested using NHTSA/DOT procedures (US National Highway Traffic Safety Administration / Department of Transportation). In this testing, the accuracy of this device was established at blood alcohol concentrations (BAC) of 0.008 and 0.032, but not at greater concentrations. Additional laboratory testing has been performed at concentrations ranging from 0.00 to 0.40% BAC.

6. Precautions

- Wait at least 20 minutes after drinking, eating or smoking before testing. Not observing this waiting period can cause inaccurate readings and damage the sensor.
- Do not blow smoke, food or liquids into the BACtrack, because this will damage the sensor.
- Do not test in areas with strong winds or smoke, or in an area where a large amount of alcohol is being consumed.
- Avoid testing in the presence of any substances that contain methyl alcohol, isopropyl alcohol or acetone. These substances may interfere with the results of the test.
- Replace the two AA batteries when the battery indicator icon reaches one bar.
- Send your tester in for periodic calibration service as required (see section 7).

- The BACtrack breathalyzer is designed to be used in a temperature range of 0–40 °C (32–104 °F).
- You cannot use the results of this product in court.
- Do not use the BACtrack as a tool to determine whether you should operate a motor vehicle or equipment, or perform any other dangerous act.
- **Do not drink and drive. Always have a designated driver when alcohol is being consumed.**

DO NOT DRINK AND DRIVE

7. Calibration

All BACtrack products are calibrated during the manufacturing process and should remain accurate for 6 to 12 months, depending on the number of tests performed and the operating conditions. If the product is used every day, it may need to be recalibrated as frequently as every month.

7.1. Calibration Instructions

Please visit <https://www.bactrack.com/calibration> or call +1 (415) 693-9756 (USA) for information on how to get your BACtrack calibrated. Outside the USA, you can also arrange calibration through your local distributor (see section 8).

8. Service and Support

If you have any questions, please contact your local store or distributor, or visit <https://www.bactrack.com>. If the breathalyzer needs to be repaired or calibrated, contact your local distributor.

BACtrack customer support (USA): support@bactrack.com, phone 877-334-6876 (Monday to Friday, 9 am to 5 pm Pacific Time).

8.1. Manufacturer

KHN Solutions, LLC

200 Green Street, Suite 200
San Francisco, CA 94111, USA
<https://www.backtrack.com>
support@backtrack.com

EU Responsible Person:

OBELIS S.A., Bd. Général Wahis 53, 1030 Brussels, Belgium
mail@obelis.net

UK Authorised Representative:

OBELIS UK LTD, Sandford Gate, East Point Business Park, Oxford OX4 6LB,
United Kingdom

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© 2025 KHN Solutions, LLC
OBELIS UK LTD, Sandford Gate, East Point Business Park, Oxford OX4 6LB, United Kingdom